**LANAFORM®**

# Stop Cell

Vibrating massager

Vibrating massage

Massage vibrant

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
Thank you for choosing the Stop Cell device from LANAFORM®. The Stop Cell device is ideal for helping you to feel better about yourself. This multifunctional massage device enables you not only to prevent unattractive fat deposits but also reduce the “orange peel” effect. Stop Cell can also be used as a body massage system.

**PLEASE READ ALL INSTRUCTIONS PRIOR TO USING YOUR «STOP CELL», IN PARTICULAR THESE BASIC SAFETY INSTRUCTIONS:**

- Only use Stop Cell for its intended use as described in this manual.
- The massage device is designed only for domestic use.
- Check that your supply voltage is suitable for the device.
- To keep the device working correctly, it is recommended to unplug it for 10 minutes after each 20 minute treatment period.
- Do not use the massage device on the following parts of the body:
  - Any part of the body which is swollen, burnt or inflamed or any part of the body where there is a rash or ulcer.
  - Aching stomach.
  - Legs affected by varicose veins.
  - Delicate genital tissue.
- People suffering from diabetes and pregnant women should consult their doctor prior to using Stop Cell.
- The massage should feel pleasant and comfortable. If you experience any pain or discomfort, please stop using the device and consult your doctor.
- This device is not designed for use by people with reduced physical, sensory or mental capabilities (including children), or people without experience or knowledge of the product, unless it may be beneficial to them and the person responsible for their safety can supervise them and has received advice and instruction in using the device. Children should be supervised to ensure they do not play with the device.
- To prevent damage to the internal components, take care not to drop the device.
- Do not get your device wet, in the bath or shower.
- To prevent overheating, do not use the device under a cover or pillow.
- Keep the cable away from hot or sharp objects.
- If the flex is damaged, it must be replaced with a suitable flex which is available from the supplier or after-sales department.

## FUNCTIONS AND USE

1. Power supply
2. Rollers
3. Charge indicator
4. 0-1 switch

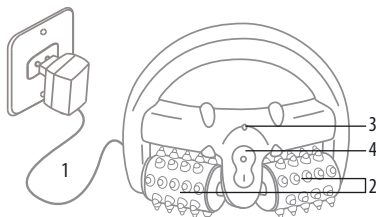Use your Stop Cell every day to prevent fat deposits and reduce the “orange peel effect” on your skin. To get the best results, combine this treatment with a balanced diet and physical exercise. The convenient Stop Cell system has been designed specifically for women who are concerned about their body.

Many women suffer from problems linked to fat deposits and water retention, which result in the appearance of cellulite. This gives their skin an “orange peel” appearance.

The “orange peel” effect is generally used to describe areas of the body where there are indentations on the skin (thighs, bums, upper arms, etc.). Sometimes cellulite is not visible to the naked eye, but only when skin is squeezed between the fingers.

## INSTRUCTIONS PRIOR TO USE

### CHARGE THE DEVICE BEFORE EACH USE

- Place the switch into the “0” position and plug the adapter into the massage device.
- The “green” indicator light will light up during charging.
- Prior to using the device for the first time, it is necessary to charge it for 24 hours. Subsequent charge time should be 7 hours.
- The operating time for the device is 25 minutes after a full charge.
- No period of use should exceed 20 minutes. Allow the device to rest for 10 minutes after each use.

### DIRECT POWER SUPPLY THROUGH THE ADAPTER

- When the battery is dead, you can continue to use the massage device by supplying power directly through the adapter provided.

- If the battery is completely dead however, you should recharge the device for 5 minutes prior to using it in this way.
- To use the system with direct power supply, simply plug the adapter from the massage device into the power supply.

## HOW TO USE YOUR MASSAGE DEVICE

This massage system has 2 soft rollers which carry out a cellulite fighting function. This system is operated through an integrated motor which generates a massage function for your muscles and skin.

For the best results, we recommend using Stop Cell regularly:


- 5 minutes each day on each area to be treated for a total of 20 minutes per session.
- The massage movement should always be carried out towards the heart.

### Notes:

The treatment may leave your skin may be slightly red and hot (this shows that blood circulation has been stimulated).

It is normal to experience these symptoms after each massage and they will disappear quickly. If you wish, you can apply a soothing lotion (not provided) after treatment.

Results may differ between women due to differences in the quantity of cellulite, skin type and lifestyle.

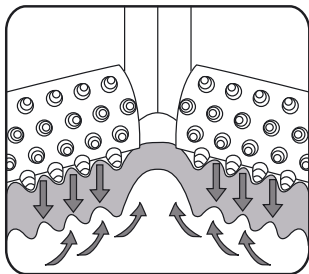Thanks to the vibrating massage, combined with the effects of skin contraction, Stop Cell really acts on cellulite.

## MAINTENANCE AND CLEANING

1. Always unplug the adapter prior to cleaning the device.
2. Never use liquids or abrasive cleaners to clean the device. Simply wipe with a damp cloth.

## ADVICE REGARDING THE DISPOSAL OF WASTE

- All the packaging is composed of materials that pose no hazard for the environment and which can be disposed of at your local sorting centre to be used as secondary raw materials. The cardboard may be disposed of in a paper recycling bin. The packaging film must be taken to your local sorting and recycling centre.
- When you have finished using the device, please dispose of it in an environmentally friendly way and in accordance with the law.

## LIMITED WARRANTY

« LANAFORM® guarantees this product against any material or manufacturing defect for a period of two years from the date of purchase, except in the circumstances described below.

The LANAFORM® guarantee does not cover damage caused as a result of normal wear to this product. In addition, the guarantee covering this LANAFORM® product does not cover damage caused by abusive or inappropriate or incorrect use, accidents, the use of unauthorized accessories, changes made to the product or any other circumstance, of whatever sort, that is outside LANAFORM®'s control.

LANAFORM® may not be held liable for any type of circumstantial, indirect or specific damage.

All implicit guarantees relating to the suitability of the product are limited to a period of two years from the initial date of purchase as long as a copy of the proof of purchase can be supplied.

On receipt, LANAFORM® will repair or replace your appliance at its discretion and will return it to you. The guarantee is only valid through LANAFORM®'s Service Centre®. Any attempt to maintain this product by a person other than LANAFORM®'s Service Centre will render this guarantee void.”

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