relax in style

MAXIWELL III

< relax in style

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< relax in style

Since being established in 2000, the name CASADA stands for high quality products on the health-spa and fitness market. Today, CASADA-equipment is available in 37 countries world-wide.

Constant further development, an unmistakeable design and a focused orientation towards a better quality of life are what characterize CASADA-products, as much as their unbeatable and reliable quality.

CONGRATULATIONS!

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal Maxiwell III.

TABLE OF CONTENTS

12	Safety instructions
12	Start up
13	Technical data
13	Equipment / Functions
15	Remote control
15	Massage techniques
16	Troubleshooting
16	Care instructions
16	Storage
16	Conformity declaration
17	Warranty conditions

SAFETY INSTRUCTIONS

Please carefully read the following instructions to before using this fitness equipment to ensure its proper function and optimal effect. Please retain these instructions for use!

- This massage device complies with the recognised principles of technology and current safety regulations.
- This equipment is maintenance free. Any necessary repairs must be performed by an authorised professional. Improper use and unauthorised repairs are prohibited for safety reasons and will void the warranty.
- Never touch the power plug with wet hands.
- Please avoid the equipment coming into contact with water, high temperatures and direct sunlight.
- Be sure to unplug the equipment from the power point during extended periods of non-use to eliminate the risk of a short circuit.
- Do not use damaged plugs, cords or loose sockets. Damaged plug or cords must be replaced by the manufacturer, service representative or qualified personnel. In the event of malfunction unplug the equipment immediately. We assume no liability for damages resulting from wrongful or improper use.
- In order to prevent excessive stress on muscles and nerves, the recommended massage time of 30 minutes should not be exceeded.
- Never use pointy or sharp objects.
- Any massage – also massage by hand – has to be avoided during pregnancy and/or in case of one or more complaints in the area of the massage: head diseases recent injuries, thrombotic illnesses, inflammations and swelling of any kind as well as cancer. Consultation of your physician is recommended prior to any treatment of afflictions and ailments.
- Do not use the device while driving.
- If you depend on electrical aids such as a pace maker, please consult your physician prior to usage in any case.

START UP

Please connect the adapter with the device and subsequently plug the power plug into the power socket. Now turn on the device by the control panel.

With the car adapter the Maxiwell III can also be used when travelling (not to be used when driving!). Simply plug the car adapter in the cigarette lighter.

TECHNICAL DATA

Measurements: 45,5 x 32,5 x 12,5 cm

Weight: 3,5 kg

Power source:

Input: AC 100-240 V ~ 50/60 Hz

Output: DC 12 V, 2500 mA

Rated power: 30 Watt

Operating time automatic: 30 W

Certificates:

EQUIPMENT / FUNCTIONS

In the new Maxiwell III the massage heads normally made from plastic are substituted with heads of genuine jade. The natural powers of jade are increased by an integrated heating system. The Massage will be intensified by sending healing warmth into the depths of the muscles. The heat is permeating through the stone slowly and is only noticeable after a few minutes. During this heating phase the tissue is already being relaxed by the rotating movement of the massage heads and surface hardening and tension are being released. The optional vibration also allows for the stimulation of the circulation and lymphatic system. Blood circulation is effectively improved and the purification facilitated.

Please position the Maxiwell III on a stable chair. The burden of your bodyweight intensifies the massage effect. Please ensure that your whole body is positioned centrally so that the massage heads run parallel to your spine.

EQUIPMENT/FUNCTIONS

The massage function will switch off automatically after 15 minutes. In order to increase the lifespan of the motors, please switch on the massage heads again only after a minimum of 30 minutes have passed. This gives the motors time to cool down again. You should also allow a period of rest for your body. We recommend conducting the massages only two, maximum three times a day in order to not overstress your muscles.

For your own safety we equip the machines with a protection against overheating. Heat is created in the motors when placed under great stress and excessive pressure. Before this could create any risk, the machines will switch off automatically and will only allow switching on again after an appropriate downtime between 30 and 40 minutes.


FOR MOBILE USE:

In addition to your Maxiwell III you will receive a 12 V adaptor for mobile use.

REMOTE CONTROL

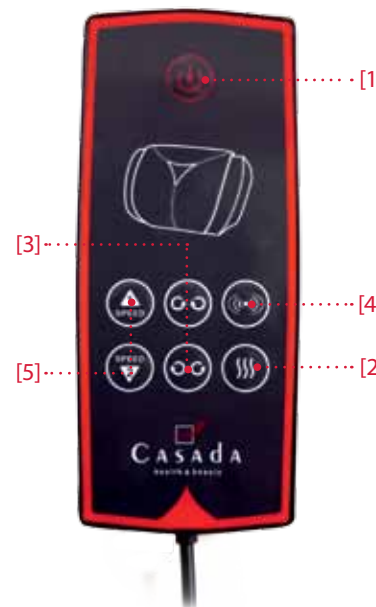- [1] On/off
Pressing the POWER, Maxiwell III stands by to start 15 minutes countdown. Pressing again or 15 minutes being up, Maxiwell III stops to work.
- [2] Heat
On/off additional heating.
- [3] Massage heads
Change direction of the rotating massage heads.
- [4] Vibration
On/off additional vibration.
- [5] Speed
Choose speed of the rotating massage heads.

MASSAGE TECHNIQUES

SHIATSU/KNEADING

Shiatsu (finger pressure) is a form of body therapy developed in Japan. The massage technique consists of gentle rhythmic twists and rotations with deep effect. Shiatsu aims to stimulate individual points in the body and to mobilise the muscles.

VIBRATION

The Vibration puts entire body regions into even vibrations, which continue all the way into the depths. Stimulating for circulation, warming and relaxing for the muscles.

HEATING

A deep and warmth-radiating effect, where one can relax perfectly. This function is especially enjoyable in case of muscle tension due to its deep effect.

TROUBLESHOOTING

Error	Adjustment
The massage device cannot operate even though it has been connected to the power supply.	<ul style="list-style-type: none">• Check the connection of the adapter.• The device could be overheated after prolonged use. It will automatically cool down to resume its normal condition. Please turn off the power and allow the device to cool off for 30 to 50 minutes.

STORAGE / CARE

Regular care will extend the life of your Maxiwell III whilst also ensuring safe operation at all times.

1. Recommended storage: dry, free of dust.
2. Keep liquids and open fire away from the Maxiwell III.
3. Clean the equipment using a dry, soft cloth.
4. Always unplug the Maxiwell III when not in use.

CONFORMITY DECLARATION

Corresponds to the European guidelines 2006/95/EC.

WARRANTY CONDITIONS

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt.

The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorized dealer. The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.

Casada International GmbH
Obermeiers Feld 3
33104 Paderborn, Germany
phone (+49) (0 52 54) 9 30 98 - 0
fax (+49) (0 52 54) 9 30 98 - 30
info@casada.com
www.casada.com

Copyright Images and texts are subject to the copyright of Casada International GmbH and may not be used for other purposes without express confirmation.

Copyright (c) 2014 Casada International GmbH. All rights reserved.

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Copyright © 2014 Casada International GmbH.