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CONGRATULATIONS

With the purchase of this massage device you have shown awareness of your health. In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal Eve.

SAFETY CONDITIONS

This equipment is maintenance free. Any necessary repairs must be performed by an authorised professional. Improper use and unauthorised repairs are prohibited for safety reasons and will void the warranty.

In order to prevent excessive stress on muscles and nerves, the recommended massage time of 15 minutes should not be exceeded.


Any massage – also massage by hand – has to be avoided during pregnancy and/or in case of one or more complaints in the area of the massage: head diseases recent injuries, thrombotic illnesses, inflammations and swelling of any kind as well as cancer. Consultation of your physician is recommended prior to any treatment of afflictions and ailments.

IN THE PACKAGE

- Massage device Eve
- Batteries 3 x 1,5 V AAA

START UP

Turn the upper cover counterclockwise, as shown in the picture, and insert three AAA dry batteries. Close the cover.

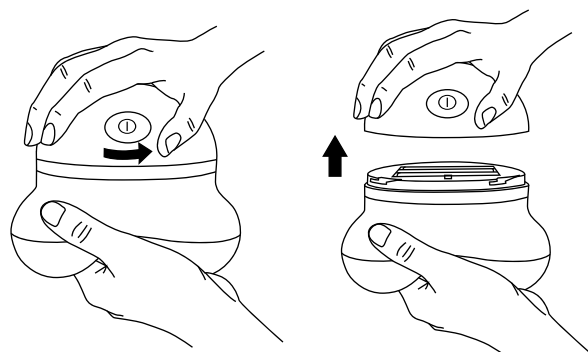TECHNICAL DATA

Measurements: 96 x 96 x 90 cm
Weight: 170 g
Rated power: 1,8 W
Battery: 3 x 1,5 V AAA
Water resistance IPX5
Certificates:

VIBRATION

EQUIPMENT / FUNCTIONS

Modern Design, a pleasant feeling and is useable everywhere. A small massage device with a big impact – that is Eve!

A precise and pleasant massage is achieved by the three illuminated massage arms and focused vibration. The device is easily and safely applied to individual areas of the body like, the back, neck, arms and legs. The massage stimulates the circulation and relieves hardened muscles and enhances blood circulation. This revitalizes the body and activates new energy.

Because of the size and the ergonomic form, Eve fits in every pocket-book and therefore can be used for example in the office, airplane and also just when travelling.

CONTROL PANEL

Press the ① button to start the massager. Press the button again to stop the massager.

TROUBLESHOOTING

Error
The massage device cannot operate.

Adjustment
Check whether the batteries are correctly installed.

Replace batteries.

The massage is too „weak“

Replace batteries.

CONFORMITY DECLARATION

Corresponds to the European guidelines 2006/95/EC.

WARRANTY CONDITIONS

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt.

The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorized dealer.

The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.

Casada International GmbH
Obermeiers Feld 3
33104 Paderborn, Germany
phone (+49) (0 52 54) 9 30 98 - 0
fax (+49) (0 52 54) 9 30 98 - 30
info@casada.com
www.casada.com

Copyright Images and texts are subject to the copyright of Casada International GmbH and may not be used for other purposes without express confirmation.

Copyright (c) 2012 Casada International GmbH. All rights reserved.