

beper

- MASSAGGIATORE ANTICELLULITE - MANUALE DI ISTRUZIONI
- ANTI-CELLULITE MASSAGER - USE INSTRUCTIONS
- MASSEUR ANTICELLULITE - MANUEL D'INSTRUCTIONS
- ANTICELLULITE-MASSAGEGERÄT - BETRIEBSANLEITUNG
- MASSAJADOR ANTICELULÍTICO - MANUAL DE INSTRUCCIONES
- SISTEM PROFESIONAL DE LIPOMODELARE - MANUAL DE INSTRUȚIUNI



Cod.: 40.500

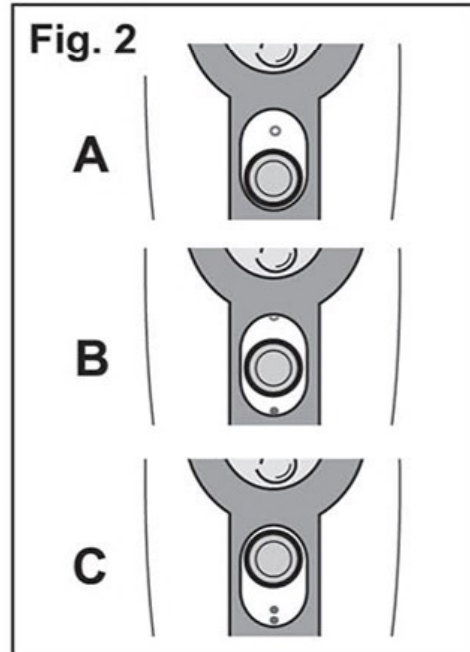
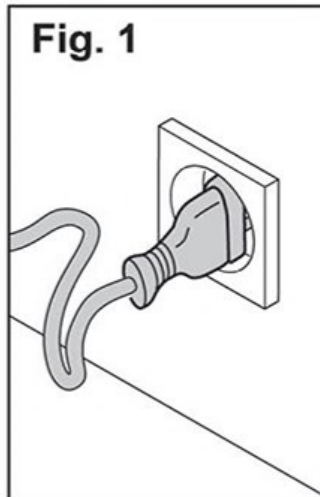
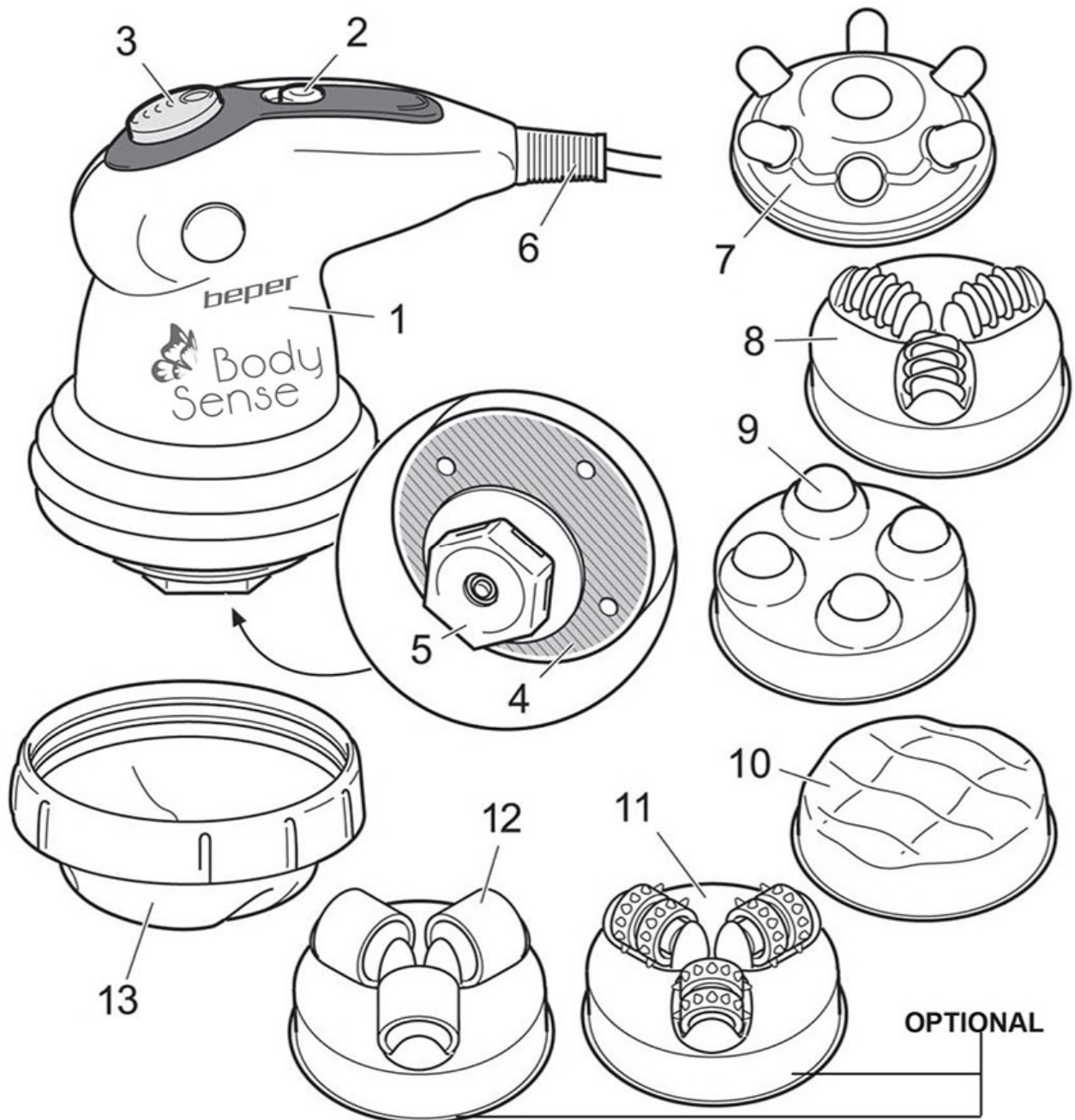


Fig. 3

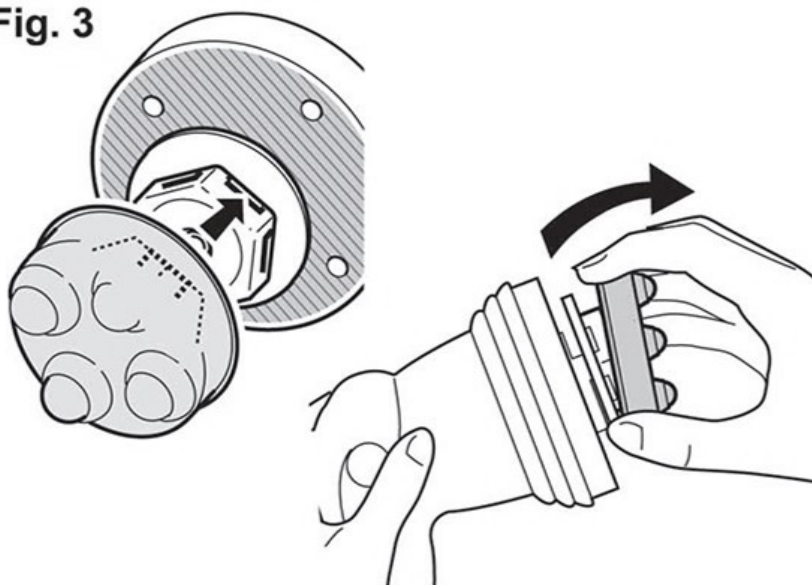


Fig. 4

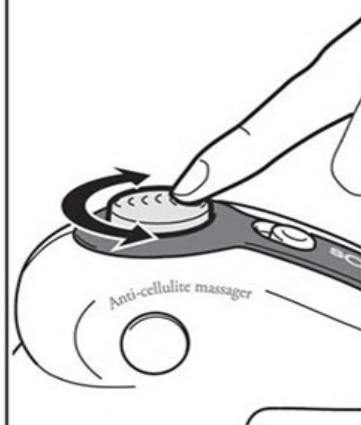


Fig. 5


5 min.



Fig. 6


5 min.


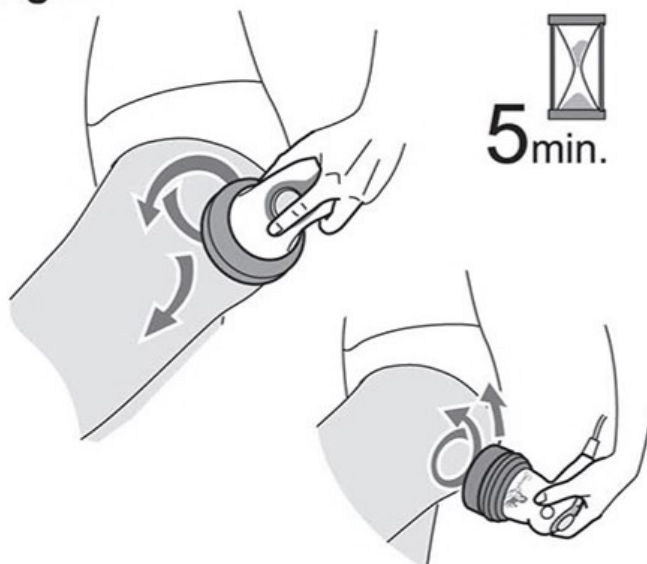


Fig. 7

2-3 min. 



Fig. 8

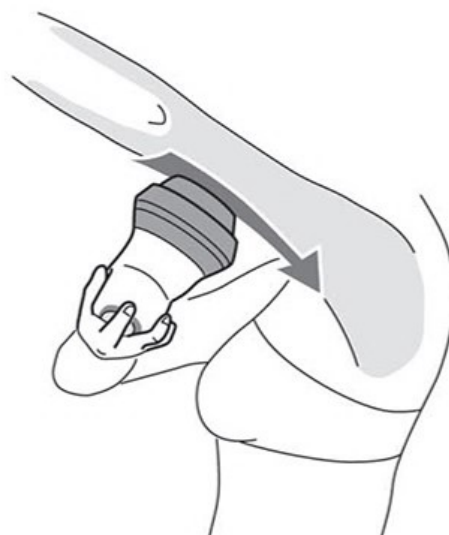

1 Safety Instructions

- Read these instructions completely before using your appliance.
- This manual is an integral part of the product. Keep it in a safe place for future reference.
- This appliance is only for domestic use, not professional use. Using it in any way other than according to the instructions could be hazardous. Do not use this appliance as a substitute for appropriate medical care.
- This appliance is not designed to be used by people (including children) whose physical, sensory or mental capacities are impaired or who lack experience or knowledge, unless they are supervised or have received instructions on how to use the appliance by a person who is responsible for their safety.
- Children must be supervised to ensure that they do not play with the appliance.
- Ask your doctor if you suffer from any health problem, for example, if you use a pace maker.
- Do not modify or repair the appliance. In there is any problem with the cord or any other part of the appliance, do not use it and take it to an authorised service centre.
- If the power cord is damaged, it must be replaced by the manufacturer or the manufacturer's after-sales service or by similarly qualified personnel to prevent possible hazards.
- Do not use parts or accessories not supplied or recommended by BEPER.
- Before plugging in the appliance to the mains, check that the indicated voltage matches the voltage in your home.
- Always unplug the appliance after use and before removing or installing parts and before performing any maintenance or cleaning task. Also unplug it in the event of a power outage. When you plug in or unplug the appliance, it must be turned off.
- Do not pull on the power cord or hang the appliance by the cord. Always unplug the appliance by pulling on the plug, never the cord.
- Make sure the cord does not come into contact with sharp edges or hot surfaces. Do not leave the cord hanging over the edge of a table or worktop, to prevent children from pulling it and causing the appliance to fall.
- Remove all elements from the inside or outside of the appliance used for protection during transport or for sales promotions, such as paper or plastic bags, plastic sheets, cardboard and stickers.

d Do not use near bathtubs, sinks, showers or other recipients containing water. If it falls into water, do NOT try to pick it up. Unplug it immediately.

- Never leave the appliance exposed to the elements.
- Do not touch the appliance with wet hands when in operation.
- Never leave the appliance working on a wet surface or cloth.

- Do not use the massager if you are suffering from any skin rash or inflammation.
- If you suffer from varicose veins, do not use the infrared light in that area, as the infrared light emits a low heat which could aggravate your circulatory problem.
- Do not use the massager in the abdomen area if you are pregnant or believe that you might be pregnant.

q **IMPORTANT:** When you want to dispose of the appliance, NEVER throw it in the rubbish bin. Take it to the CLEAN POINT or waste collection centre closest to your home for processing. You will thus be helping to care for the environment.

2 Description

Main parts

- | | | | |
|---|------------------------------------|----|--|
| 1 | Main body of the appliance | 8 | Spiky roller massage head |
| 2 | Off/massage/massage + IR switch | 9 | Ball massage head |
| 3 | Variable massage intensity control | 10 | Wavy surface massage head |
| 4 | Infrared light emitter | 11 | Parallel spiky disk massage head |
| 5 | Massage head adaptor | 12 | Flat roller massage head |
| 6 | Cord | 13 | Protective fabric for wavy surfaced massage head |
| 7 | Finger massage head | | |

Electrical characteristics

Charger power supply voltage: 230-240V, 50Hz, 25W

Professional Lipomodelling system.

Fights cellulite and flaccidity efficiently at home

This latest-generation massager allows you to enjoy an effective therapeutic massage at home and is especially recommended for people suffering from cellulite, fluid retention and flaccidity. Continuous use of the apparatus will help you to:

- Improve blood circulation.
- Obtain firmer and smoother skin.
- Eliminate fluids and toxins.
- Eliminate accumulated fat and "orange peel" skin.

3 Selecting the Massage Head

Head Types and Functions

- Finger massage" massage head (7) and "Flat roller" massage head (12)
For use on painful cellulite and muscular pains.

! Do not use on the back.

- "Spiky roller" massage head (8) and "Parallel spiky disk" massage head (11)

This head can be used to eliminate cellulite, promote firmness, stimulate blood flow and return, tone muscles, drain the lymphatic system, reduce size, warm up the muscles, provide relaxation, relieve stress, etc.

! Do not use on the back.

- “Massage balls” head (9)

Stimulates and tones the muscles, helping to firm the skin. It improves circulation in the abdominal area, eliminating localised toxins and improving the intestinal tract, preventing constipation.

- “Wavy surface” massage head (10)

Tones and improves blood circulation. For use in the thigh and buttock areas. For greater comfort, cover the wavy surface massage head (10) with the protective fabric (13).

- Use of the “head protector” cloth head (13)

- ☞ The cloth accessory (13) can only be used with the massage head (10).










“Wavy surface”









Types of Cellulite

Identify the type of cellulite you suffer from, to help you choose the most effective massage head.

- Oedematous Cellulite: Feels spongy, it does not hurt and moves when touched or when you change position.
- Painful Cellulite: Pain is felt when the area is pressed.
- Cystic Cellulite: Grainy and hard to the touch. Usually accompanied by stretch marks. Its temperature feels lower compared to the rest of the body.

Treatment Type

ZONE	CONDITION	HEAD	INTENSITY
LEGS AND BUTTOCKS	Oedematous Cellulite Flaccidity Cystic Cellulite	1st STEP 	
		2nd STEP 	
		3rd STEP 	
	Painful Cellulite Tiredness		
			
	Fluid Retention		

ZONE	CONDITION	HEAD	INTENSITY
ABDOMEN	Oedematous Cellulite Painful Cellulite Flaccidity Constipation		
	Cystic Cellulite		
ARMS	Oedematous Cellulite Painful Cellulite Flaccidity		
	Cystic Cellulite		






4 General Recommendations for Use


- The intensity control allows you to regulate the intensity of the massage as required. Start at a gentle speed and progressively increase it, depending on your sensitivity to the massage.
- If using special treatment creams (firming and anti-cellulite creams, for example), apply these after the massage, which prepares the skin and its tissues to absorb the ingredients.
We always recommend you use a moisturising lotion after the massage. You can apply the cream using either the "wavy surface" head (10) or the "flat roller" head (12).
- Do not apply creams when using the protective fabric (13) on the "wave-shape" head.

5 Recommended Treatment Duration

- Apply a 2-3 minute massage on the abdomen and a 5 minute massage on the legs and buttocks. You may establish the duration of the massage depending on your sensitivity. The total treatment must be approximately 15 minutes.
- Do not use the appliance for more than 15 minutes. Using the appliance continuously for more than 15 minutes may cause it to overheat, reducing its life. Let it cool down before continuing use.
- Use Sculptural every day until you notice the results, which may be subsequently maintained by using it three days a week. To maintain the results over the long-term, adapt your lifestyle habits and food.

6 Use

-  This appliance has been designed to provide therapeutic massages to help prevent cellulite, fluid retention and flaccidity. Do not use it for any other purpose.
-  Massaging with this professional appliance generates a thermal effect in the treated area; this thermoactive effect helps to mobilise the accumulated fats and toxins. To gently intensify the treatment, the appliance has an infrared light that slightly increases the thermal action of the appliance.
 - 1 Plug in the appliance (fig. 1).
 - 2 Move the switch (2) to the massage position (fig.2 B) or to the "massage + infrared position (IR) position" (fig.2 C) according to the desired treatment.
 - 3 Fit the correct head, (7), (8), (9), (10), (11) or (12), for the treatment you require to the adaptor (5) (fig. 3).
-  To remove the head, pull firmly at the edge to release the fastening tabs. (fig. 3).
-  To choose which massage head to use, follow the instructions in the "Treatment Type" table.
-  You may alternate the use of different heads during the same massage.
 - 4 Set the intensity of the massage using the intensity control (3) and place the massager on the area on which you are going to work.
 - 5 Set the intensity of the massage using the control (3), according to the area to be treated and its sensitivity to massage (see the "Treatment Type" table).

To do so, turn the control to one side or another with your finger, depending on the desired intensity (fig. 4).
 - 6 Hold the massager firmly with both hands and move it over the skin in the direction of the blood circulation (see "Massage Direction").
-  Without pressing too hard, simply glide the massager over the area to be treated.
- 7 When you have finished using the massager, move the switch (2) to the "off" position (fig.2 A).
- 8 Unplug the appliance.

7 Massage Direction

Always glide the massager over the skin in the direction of the blood circulation.

- Legs

Place the appliance on the thigh, close to the knee, and glide it gently upwards towards the hip (fig. 5).

- Buttocks

Combine upwards and circular movements. You can end an upward massage on the legs with downwards and circular movements on the buttocks (fig. 6).

- Abdomen

Place the appliance on the navel and move it gently outwards (fig. 7).

✋ Never make circular movements around the navel.

- Arms

Place the appliance close to the elbow and glide it gently towards the armpit. Repeat the massage from the armpit towards the elbow (fig. 8).

✋ NEVER apply a massage on a bone (knee, ankle, elbow). The protruding parts of the massage head could knock the bone, causing an injury.

✋ Only use the massage heads indicated in the "Treatment Type" table.

8 Cleaning and upkeep

✋ Do not use solvents, detergents or abrasive products.

1 Unplug the appliance.

2 Use a dry cloth to clean the outside.

Storage

- When the appliance is not being used, unplug it and store it in a cool, dry place.

- To prevent the cable from becoming damaged, do not twist it, pull it or wind it around the appliance. If it is damaged, have it replaced immediately by qualified personnel.

TECHNICAL DATA

Power 25W

Voltage 230V ~ 50Hz

For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.



The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment. The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

GUARANTEE CONDITIONS

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product

TECHNICAL ASSISTANCE

Even after the guarantee period we will always pay attention to repairs for defected appliances.

For technical assistance and/or repairs after the guarantee period you can directly contact the below address:

CONTACT YOUR DISTRIBUTOR IN YOUR COUNTRY OR ATER SALES DEPARTMENT BEPER. E-MAIL assistenza@beper.com WHICH WILL FORWARD YOUR ENQUIRIES TO YOUR DISTRIBUTOR.