


DISPLAY




DISPLAY View baby's weight and other information



 **ON/OFF button** Switch ON or OFF the scale.

[TARE] button Allows you to deduct initial tare weight.

[UNIT] button Allows you to change weight unit.

 **button** Allows you to use AUTO HOLD and MANUAL HOLD function, explained in this manual.

NORMAL WEIGHING FUNCTION


1. Press " button to switch ON the scale.
2. "0.000 Kg" and battery icon will show. "0.00 lb" will show if lb unit is selected.
3. Put the baby onto the scale. The weight is determined when " appears.

NOTE: Since a baby is a moving object, the weight on the display will fluctuate during weighing process. It will take about 15 seconds for the scale to catch a weight of the moving baby and lock it. The locked weight will show and then the scale will switch off automatically).

4. "EEEE" will show if the scale is overloaded.

TARE FUNCTION

This feature allows you to weigh and to exclude the weight of any tissue or blanket placed on the scale when weighing the baby naked.

1. Press " button to switch ON the scale.
2. "0.000 Kg" and battery icon will show. "0.00 lb" will show if lb unit is selected.

SPECIFICATIONS

Range from 0,02kg to 18kg / from 1oz to 40lbs

Graduation 5 g/oz

Power supply 2 x batteries AAA (3.0 VDC)

Battery indicator Battery icon will show "Lo" if battery power is down.

Display Backlight LCD panel

Display dimensions 74 x 28,5mm

Dimensions 560 x 330 x 46mm

Net weight 1.68kg



Storage environment -10°C – 55°C

Operating environment 10°C – 40°C


Precision of weight measurements 0.005g / 0.01lb

BATTERIES INSTALLATION

1. Turn the scale and remove the battery cover which is located at the bottom.
2. Insert two standard batteries size AAA being sure to observe the polarity indicated in the battery compartment.
3. Close the batteries cover and place the scale on a flat and hard surface to weight the baby. Carpet or soft surface is not suitable for weighing.
4. Before putting a baby on the scale, be sure that the surface of the scale is completely dry and not slippery.

3. Place the blanket on the scale. Weight of the blanket will show.
4. When “” appears press the [TARE] button immediately.
5. “0.000 Kg” (0.00 lb will show if this unit is selected) and TARE will show.
NOTE: Leave the blanket on the scale. If you remove the blanket from the scale, the display will show a negative weight.
6. Put the naked baby onto the blanket. Net weight of the naked baby will show when “” appears. The scale will switch off automatically after 3 minutes.

UNIT SELECTOR FUNCTION

1. Press “” button to switch ON the scale.
2. Before putting the baby onto the scale, press [UNIT] button to select the weight unit.
3. Determined weight unit will show when the scale is switched ON next time.

WARNING AND GENERAL SAFETY RULES

- Read the instructions of this manual before using this appliance. It contains important information regarding the use of this appliance and safety rules.
- Any incorrect operation and improper use may damage the appliance and cause danger to the user.
- Never leave the baby unattended during weighing
- Use the scale only to weight babies below the maximum capacity. Over loading will cause damage and danger.
- To clean the balance use a soft cloth dampened with a bit of water. Do not immerse the scale in liquid or rinse under running water.
- For cleaning do not use detergents or solvents aggressive or abrasive
- When the scale is in use, keep away from any electromagnetic wave source such as microwave oven, mobile phones and wireless network.
- Do not short circuit the battery poles.
- Do not try to repair the appliance by yourself to avoid any dangers related to electricity.
- Never let children play with any packing materials, especially plastic bags. Plastic bags may cause suffocation.
- In case of an electrostatic malfunction, reinstall batteries to resume normal performances.
- If the display shows the error message “Err”, switch OFF the scale. Then switch ON again and follow the instructions.
- Never leave anything onto the platform of scale when not in use.
- Do not use the balance to a sustained load as it can damage the mechanism of measurement.
- Stored weight record will be deleted when new batteries are inserted.
- The scale will switch OFF automatically in 1 minute after the scale is switched ON and not in use. To switch OFF the scale immediately, press “ ” button and hold about half second.
- Avoid abrupt movement or that the scale falls.
- The scale is intended for private use and to measure the weight of the baby. It is not allowed to use the scale in medicine or for purposes other than private use.
- As the scale is intended exclusively for private use is exempt from the requirement to submit the instrument to the periodic inspections of the law, stating that they are to protect the accuracy and reliability of the regularity of transactions.



Dear customer,

Thank you for choosing this baby scale. Please read these instructions for use carefully and keep them for future use.

SPECIFICATIONS


- Accurately measures weight in 5g increments up to 20kg.
- HOLD feature captures precise weight measurement of an active baby.
- TARE feature subtracts added weight of diaper, clothing or blanket.
- Easy-to-read, white backlit LCD display.
- Unique touch button design.
- Includes 2xAAA batteries.
- Auto power off and low battery indication.

SPECIFICATIONS




POWER SUPPLY

Please insert batteries or remove the insulating foil from the batteries before using. When "Lo" appears in the display, the batteries need to be replaced.

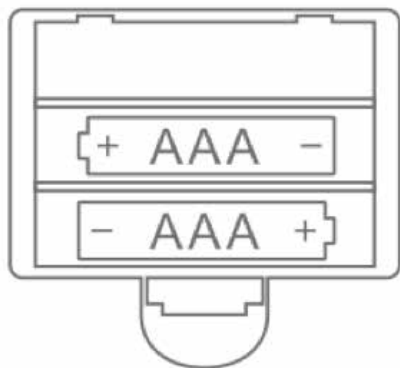

POWER SUPPLY


The weight is displayed in steps of 5g from a minimum weight of 30g. Press the button  to switch on the scale.

When "0.000kg" appears on the display, the scale is ready for measurement. Place the baby on the scale carefully. Take particular care with the baby's head.





Hold Function:

The scale has a **HOLD** function: **AUTO HOLD**, **MECHANICAL HOLD**.

AUTO HOLD:

When the baby is placed on platform, they may move around. When they are still enough, the AUTO HOLD function will establish the weight on display. The sign  flashes and then weight is established.

MECHANICAL HOLD:

If the baby is moving around too vigorously, you may select the MECHANICAL HOLD function by pressing  button. The sign  flashes and weight is established on display. The sign  remains on display until you press  button or "TARE" button. Then it shows "0.000kg" which means the scale is ready for next weighing.

Tare Function:


If you prefer to place a towel or blanket on the scale and weigh the baby, please do so BEFORE switching the scale on. If the towel or blanket weighs over 30g, please press TARE button.

"TARE" is shown on display. Then place the baby onto towel or blanket. It shows baby's weight.

Units:

The scale is with units of kilograms (kg) and pounds (lb), which can be switched by pressing "UNIT" button.

SWITCH OFF

After weighing is finished, you can switch off the scale by long pressing  button. Or, the scale will automatically switch off if no weighing is carried out within 180 seconds.

OVERLOAD INDICATION

The scale can accept a maximum load of 20 kg. When overloaded, “EEEE” sign will appear on display.

CARE AND ATTENTION

- Place the scale on a hard, flat surface to ensure greatest accuracy.
 - Cleaning: You can clean the scale with a damp cloth and a little washing up liquid if necessary. Never immerse the scale in water or rinse it under running water.
 - Protect the scale against knocks, damp, dust, chemicals, drastic changes in temperature, electromagnetic field and source of heat that are too near.
 - Keep packaging material away from children.
 - DO NOT dispose of the old batteries in fire. The batteries may explode or leak. Recycle or manage them as hazardous waste.
 - If the scale is not used in a long period, please take the batteries out.
-

WARRANTY - TERMS AND CONDITIONS

This product benefits of a warranty of 24 months on material and manufacturing defects, starting from the date of purchase (see sales slip).

The warranty of 24 months does not include damages caused by the usual use of parties identified as “consumable” (e.g., batteries, brush heads or parts subject to usury).

The legal guarantee of 24 months is void if:

1. The product has undergone aesthetic damage due to improper use not in accordance with instructions in the manual.
2. This product has been modified and/or tampered with.
3. The cause of the failure was due to poor maintenance of the individual components and/or accessories and/or supplies (e.g. oxidation and/or scaling due to the retention of water or other liquids, sediment blocking the sensor, leak of corrosive liquid from batteries).

The following is excluded from the legal guarantee of 24 months:

1. Costs related to replacements and/or repair of parts subject to wear or costs for ordinary maintenance of the product.
2. The costs and risks involved in transporting the product to and from the store where you purchased or otherwise authorized collection centre to receive the products under warranty.

3. Damage caused by or resulting from improper installation or improper use not in accordance with the directions in the instruction manual.

4. Damage due to natural disasters, accidental events or adverse conditions not compatible with the product.

5. Defects that have a negligible effect on product performance.

The manufacturer, distributor and all the parties involved in the sale do not assume any liability for losses and economic damage from any malfunction of the product. In accordance with current regulations the manufacturer, distributor and all the parties involved in the sale are not responding in any case for damages, including direct, indirect ones, loss of net income, loss of savings and additional damage and other details consequences going beyond the damage caused by the breach of warranty, contract, strict liability, wrongdoing or due to other causes, resulting from the use or inability to use the product and/or paper and electronic documents, including the lack of service.

For further information on the help service visit the website **www.nuvitababy.com**.